

LAARMANN[®]

Innovators in Solids



LM-CM300 CUTTING MILL USER MANUAL

Contents

1	Safety information	4
1.1	Safety symbols	4
1.2	Safety instructions	4
1.3	General safety Instructions	6
2	General information	7
2.1	Repairs and modifications	7
3	Technical features	8
3.1	Construction	8
3.2	Technical data	8
3.3	Manufacturer details	8
3.4	Contact details	8
4	Transport	9
5	Overview of cutting mill	10
6	Installation	11
6.1	Unpacking	11
6.2	Selecting the correct place	11
6.3	Connecting the power cord	11
6.4	Environment conditions	11
7	Working instructions	12
7.1	General working instructions	12
7.2	Starting the process	12
7.3	Stopping the process	12
7.4	Residual risks	13
7.5	Intended use	13
7.6	Inproper use	13
8	Workingstations	14
8.1	Installation Area	14

Contents

8.2	Safety and Environmental Requirements	14
8.3	Working Environment	14
9	Operating the control panel	15
9.1	Info screen	15
9.2	Overview and status	15
9.3	Settings menu	15
9.4	Timed mode and revolution mode	16
9.5	Service menu	17
9.6	Errors and notifications	17
10	Operating the CM300	18
10.1	Insert the bottom sieve	18
10.2	grinding sample material	19
10.3	Replacement static knives	19
10.4	Adjust the gap	20
10.5	Adjusting by setting the rotor knife	20
11	Preparation before grinding operation	21
11.1	Standard hopper	21
11.2	Start the instrument	21
11.3	Stop the instrument	21
12	Troubleshoot and Maintenance	22
12.1	Troubleshoot	22
12.2	Maintenance and lubrication	23
12.3	Wear	24
12.4	Cleaning	24
13	CE certificate	25

1 Safety information

Before using the machine, make sure to read and understand this manual thoroughly. Keep the manual close to the machine, easily accessible to all the users. Improper operation can cause injury to persons or damage to the equipment.

1.1 Safety symbols

The following are the warning symbols that are used in this manual. Failure to comply with the sections of text marked in this way result in danger. Injuries are possible. Damage to the equipment and accessories is possible, special care is required.



This symbol indicates a potential risk and alerts you to proceed with caution.



This symbol indicates the presence of high voltage and warns the user to proceed with caution.



This symbol indicates risks associated with hot surfaces.

1.2 Safety instructions



Safety instructions

Claims for damages in any form whatsoever, for injury to persons or damage to the machine, caused through non observance of the following safety instructions, are excluded.



Use according to the intended purpose

Do not make any alterations to the machine and use only approved spare parts and accessories. Otherwise the Declaration of Conformity will lose its validity and this will also lead to the loss of any guarantee claims.



Transport

Do not knock, shake the machine during transport. Otherwise the electronic and mechanical components may be damaged.



Packing material

Please keep the packing material for the duration of the guarantee period. In case of a complaint and return of the machine in unsuitable packing material, your guarantee claim will be lost.

**Temperature variations**

If the CM300 is subjected to high temperature variations, protect it against condensed water. Otherwise the electronic components may be damaged.

**Ambient temperature**

If the temperature drops below +5°C or exceeds +40°C, electronic and mechanical components can be damaged. Performance can be changed to an unknown extent.

**Atmospheric humidity**

If the humidity exceeds 85%, electronic and mechanical components can be damaged. Performance can be changed to an unknown extent.

**Electrical connection**

If the values for the mains power supply on the name plate are not observed, the electrical and mechanical components may be damaged.

**Materials**

Observe the relevant regulations and directives for handling chemicals and hazardous materials. grinding or milling of materials, which give a risk of fire or explosion, is prohibited.

**Cleaning**

Do not clean the CM300 under running water. Danger to life through electric shock. Use only a soft cloth moistened with water. Cleaning agents and solvents should not be used, not for cleaning the dividing tools either.

**Repair**

For your own safety, repairs must be carried out only by authorized service technicians.

1.3 General safety Instructions

- Only instructed and authorized personnel should use this equipment.
- Wear protective clothing when operating or cleaning the machine. This includes safety glasses, safety shoes, ear protection, and a dust mask.
- When working above 2m and not on a platform use a safety harness.
- The equipment runs in an continuous process, during this only the process operator and authorized helpers should be near the equipment.
- The equipment is supplied with lethal voltages, only authorized personnel must work on or open any electrical equipment on the CM300.
- Only use correct tools to adjust any of parts of the equipment, if these are not supplied then consult with LAARMANN as to what can be used instead.
- Only use spare parts and accessories approved by LAARMANN
- Only use the specified grease or lubrication method and lubricant, when lubricating the machine.

2

General information

The CM300 can process a large varieties of materials and at the end of each milling, the door can be opened, the rotor with cutting knives can be removed by hand very easily to clean up of the grinding chamber.

The CM300 can be used for batch and continuous grinding of samples

The CM300 can mill flexible, medium-hard and fiber products and mixed products.

The unique shape of the milling tool drive with variable speed can be exempted from the effects of grinding material for fast, efficient grinding. Co-products. The milling tools do not apply to grind damp or wet materials.

The package contains:

- LM-CM300 Cutting Mill (including the controller and funnel)
- Rotor, knife set and bottom sieve
- Collection bucket: 1
- Power cable: 1
- Key to access the electrical housing
- Electrical drawings (placed inside control box or housing)

2.1

Repairs and modifications

This manual does not contain any repair and modification instructions. For safety reasons, repairs and modifications may only be carried out by LAARMANN authorised personel or by qualified service technicians.

When repairs are required reach out to LAARMANN.

3 Technical features

3.1 Construction

The housing of the CM300 is generally made of steel with powder coated RAL5005 paint.

The CM300 is equipped with a safety cover and closed down with a switch. Starting the machine is possible only when:

- power is connected;
- the emergency switch is not activated;
- all parameters are given to the digital program of the machine.
- the door is closed

3.2 Technical data

Volume collector	5 litre
Rotation speed	750-3018 rpm
Power supply	230V 50Hz
Ambient temperature	5°C – 40°C
Atmospheric humidity	< 85% RH
Dimensions (w x l x h)	600 x 684 x 925 mm
Weight	90kg

3.3 Manufacturer details

Manufacturer	LAARMANN Group B.V.
Model	CM300 Cutting Mill
Serial Number	
Standard	CE

3.4 Contact details


LAARMANN Group B.V.	+31 6 23 40 00 33
Op het Schoor 6	www.laarmann.eu
6041 AV Roermond	info@laarmann.eu
Netherlands	

4

Transport

On delivery, the device is screwed onto the transport pallet. Remove the screws from the transport pallet. Using a forklift, lift the CM300 up. Always use the forklift to lift and move the CM300. Placing the CM300 close to the installation point, prior to unpacking, will prevent unnecessary moving and therefore reduce the risk of injury. In summary, please follow the next steps after receiving the machine

- Ensure a visual inspection upon delivery and take note of any damage that may have taken place.
- Remove the screws that hold the front of the packing crate in place, release the front of the crate and set aside out of the way.
- Remove the top, sides and rear panels of the packing crate.
- Carefully remove the plastic protective wrapping and discard.
- Carefully cut the strapping used to secure the CM300 to the crate base and discard.

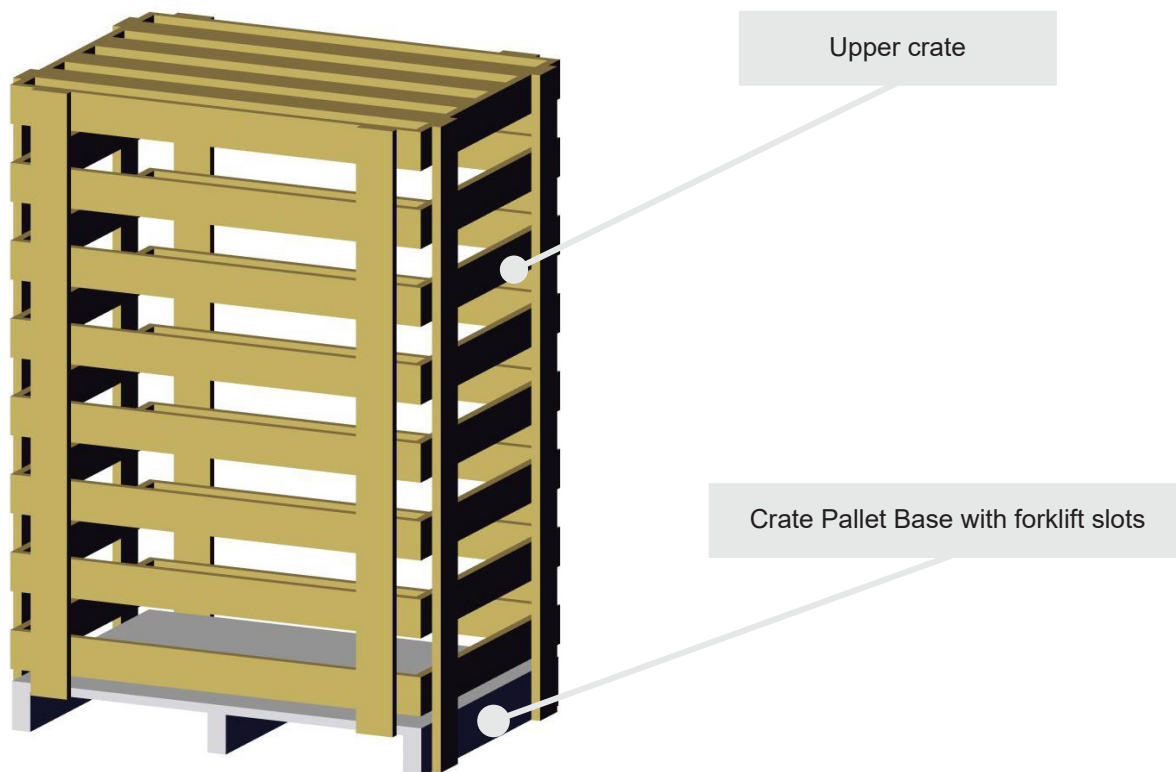

5

Overview of cutting mill



Element	Description	Function
A	Hopper	Used to feed the material safely into the cutting chamber
B	Rotor knives	Cutting the sample by rotation
C	Stator knives	Cutting the sample
D	Bottom sieve	Control the particle size at the bottom sieve plate material
E	Grinding chamber door	Access to the grinding chamber
F	Collection barrel	Receiver samples
G	Operating panel	Controlling equipment
H	Control cabinet	Installation of electrical control unit
I	Enclosure handle	For opening and closing the door of the grinding chamber
J	Support frame	Installation of the bare unit
K	Power switch / socket	Power supply to the instrument
L	Back plate Control cabinet	Closure control cabinet

6 Installation

6.1 Unpacking

Before the installation, carefully examine the delivery for possible damage or missing parts. When applicable open the box and take the equipment out of the box.

Check that the machine has not been visibly damaged during the transport. Please keep the packing material for the duration of the guarantee period. In case of a complaint and return of the machine in unsuitable packing material, your guarantee claim will be lost.

If any kind of damage occurred during transport, immediately make a complaint to the carrier. Any incorrect delivery or missing parts should be reported to the supplier.

6.2 Selecting the correct place

When selecting the right place for the machine, please consider the following:

- Put the device on smooth, stable and level concrete floor.
- Leave enough space around the device, that you will easy control and maintain it.
- Don't use the device in surroundings, where there are fast temperature and humidity changes.
- Avoid places exposed to direct sunlight and places nearby heating devices.
- Avoid places, where the possibility of shocks and vibrations exists.
- Bolt the frame securely to a floor that can safely support the weight of the unit using anchors if applicable
- Place the machine within 3 meters of a electricity socket.

6.3 Connecting the power cord

The correct voltage and frequency for the CM300 are given on the name plate. Ensure that these values correspond to the available power supply system. Connect the power cord to a grounded wall socket.

To avoid interference from noise, surges and spikes, a dedicated line is preferred. If no such line is available, avoid lines to which powerful electric motors, refrigerators and similar devices are connected.

6.4 Environment conditions

The machine has been built for operating in process and laboratory environment.

Therefore the environmental conditions should be the following:

- Temperature from 5°C to 40°C
- Humidity up to 85% RH, non-condensing

7 Working instructions

7.1 General working instructions

The CM300 is a high performance product. Because of the robust execution, the CM300 is a machine with many different application possibilities in laboratories.

7.2 Starting the process



The grinding process can only be started, when the door is closed and the power is connected



Do not open the doors during the dividing process, also not by force.



Do not mill inflammable or explosive samples! For these kind of samples, please contact LAARMANN for special features.

7.3 Stopping the process

The following actions can stop the process.

- Pressing the stop button on the display
- In case of emergency the emergency stop button can be pressed
- The cutting process automatically ends when the grinding time elapses. The process can only be restarted after the grinding time has been reset.

7.4 Residual risks

Despite correct installation and operation, certain residual risks cannot be fully eliminated. These include:

- Mechanical hazards: Risk of injury from moving parts if the safety cover is not properly closed or interlocks are bypassed.
- Noise exposure: Prolonged operation can generate sound levels that may require hearing protection.
- Dust or aerosol release: Depending on the material, fine particles may escape if no dust extraction system is installed.
- Hot surfaces: After extended operation, the rotor and knives housing may become hot.
- Handling risks: Lifting and loading of the cutting mill cause strain or injury if performed without appropriate lifting aids.
- Electrical hazards: Possible risk of electric shock in the event of damaged cables or improper grounding.
- Operators must always observe the safety instructions, wear suitable personal protective equipment (PPE), and follow the maintenance and cleaning procedures to minimize these risks.

7.5 Intended use

The cutting mill is designed for the size reduction of soft to medium-hard, fibrous, and tough materials through a cutting and shearing action. It is suitable for processing dry or slightly moist samples into uniform particle sizes for subsequent analysis, recycling, or further processing.

- Typical materials include:
- Plastics and polymers
- Rubber and elastomers
- Textiles and fibers
- Wood, paper, and cardboard
- Plant materials and agricultural products
- Leather and similar organic materials

The device is intended for use in laboratories, pilot plants, and industrial environments by trained personnel.

7.6 Improper use

Any use other than that described in this manual is considered improper and may lead to personal injury or damage to the machine. Improper use includes:

- Processing explosive, highly flammable, or reactive substances
- Milling materials harder than the specified limits
- Operating the machine without installed safety devices
- Using the mill for wet grinding unless explicitly approved

Use of the cutting mill for purposes other than those described may result in equipment damage, safety hazards, and voiding of warranty.

8 Workingstations

8.1 Installation Area

The cutting mill must be installed in a well-ventilated, clean, and dry area that meets the following conditions:

- Stable and level floor capable of supporting the total weight of the machine
- Sufficient space around the machine (minimum 1 meter clearance on all sides) to allow safe operation, maintenance, and knife/rotor replacement.
- Non-slip flooring to prevent accidents during loading and unloading.
- No exposure to vibration sources or excessive dust from other equipment.

8.2 Power Supply and Utilities

- The workplace must have the correct electrical supply according to the machine's rating plate (see Chapter 3.2 – Technical Data).
- Provide easy access to the main power switch for emergency disconnection.
- If the machine includes a cooling or exhaust system, ensure proper connection to water, air, or ventilation ducts as specified.

8.3 Safety and Environmental Requirements

- The area must be free from flammable materials and adequately lit.
- Operators must have access to personal protective equipment (PPE) including safety glasses, hearing protection, gloves, and protective clothing.
- Emergency stop buttons and first aid equipment must be easily accessible.
- If grinding produces dust or fine particles, a local exhaust or dust extraction system must be installed to maintain safe air quality.



8.5 Working Environment

- Only trained personnel are allowed to operate the cutting mill.
- The workplace should allow safe handling of the knife sets and rotors
- Adequate space must be provided for sample handling, weighing, and cleaning operations.

9 Operating the control panel


9.1 Info screen



When pressing the logo in the top left you will be redirected to the information screen which includes contact details and machine title.

For changing the language you can contact us on info@laarmann.eu

9.2 Overview and status

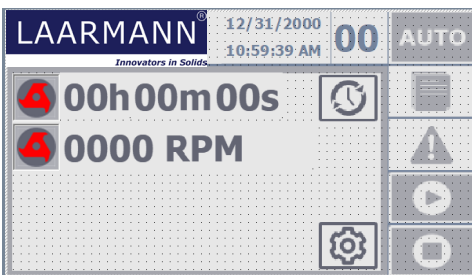

This is the overview screen of the control panel. From here, all the options can be accessed and the machine can be operated. The remaining time, current rpm and the revolutions are displayed here; depending on the mode you are using

The border color around the interface indicates the status of the machine. Green means the machine is on, red means there is an error and grey means the machine is off



Press this icon to reset the timer



Press this icon to view the settings of the program



Press the exclamation mark icon to navigate to the notification page to review current errors and notifications.



Press this icon to start the program



Press this icon to stop the program

9.3 Settings menu




In the settings menu the running time, amount of RPM (revolutions per minute) and amount of time can be adjusted.

Press on the number on the screen and select the amount. Press the enter key to confirm.

9.4


Timed mode

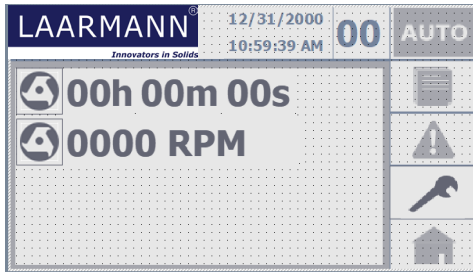

Time Mode

In Time Mode, you can set the running time required for a complete operation.

When a desired time is selected, the revolution counter (REV) is automatically set to 0, and the program will operate in Timed Mode once the machine is started from the Main Screen.

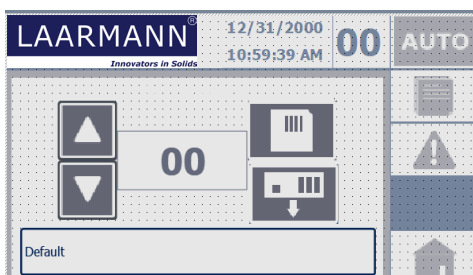

Press this button to view the recipe selection page

This is an overview of the recipe selection page for automatic mode. In this page all the settings that you have entered can be saved or loaded. The initial indicated number corresponds to the current active recipe/program.



Load current selected program into the machine.



Save active settings in the selected program number.






Use the arrows button to select other programs



Press this button to return to the main page

9.5

Service menu

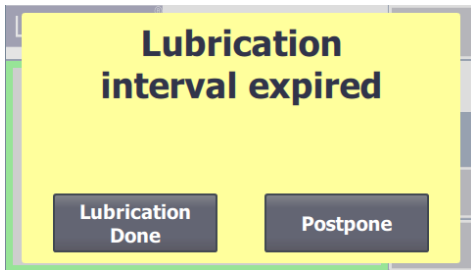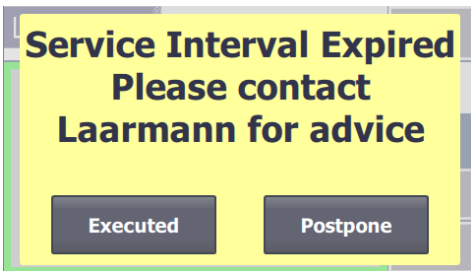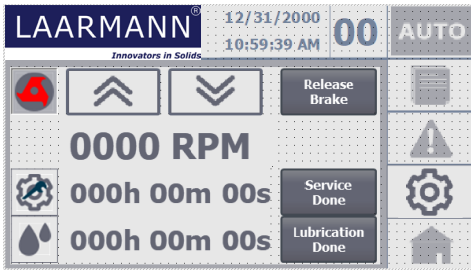

When the lubrication and/or interval or service interval ends, the corresponding warning screen appears. The machine can continue to operate without interruption. It is only a warning.

The operator chooses whether he wants to perform the lubrication/service at that moment, or whether he postpones it. With "postpone" the interval will be halved each time, until this message will eventually appear every 10 minutes. Performing or postponing is registered in a log file.

Once the postponed lubrication or service has been completed, confirm it via the maintenance screen.


Pressing Lubrication/Service done will reset the timer to 120/72 hour(s)



Use the arrows button to slowly rotate the drive of the cutting mill form maintenance

9.6

Errors and notifications

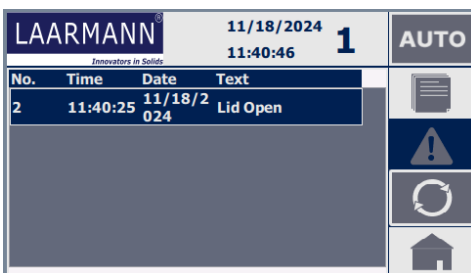

On the notification page there can be a list of errors that need to be solved before continue using the machine. Press the error row to view the full error description



Press the reset button to reset all fault messages. If there are any unresolved faults, its corresponding faultmessage will re-appear.

Warnings (i.e. service/lubrication message) can not be reset by this button, but then again they will never interlock the operation of the machine.

10

Operating the CM300

- Only when the CM300 is stopped the door can be opened
- Pressing the STOP button on the control panel, the instrument stops running
- Press the **red button (R)** and pull the lever with the **red handle (S)** to open the access door
- WAIT WITH OPENING THE DOOR WHEN MACHINE HAS COMPLETELY STOPPED



S

R





When the motor is running do not open the door of the grinding chamber! There is risk of personal injury!

10.1

Insert the bottom sieve

- Turn off Power switch CM300
- Open the door grinding chamber
- Insert the bottom sieve into the grinding chamber

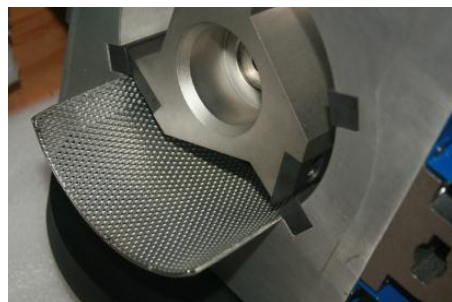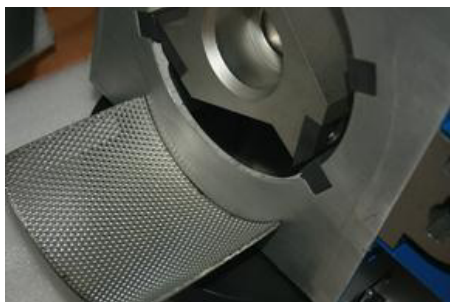




10.2 grinding sample material

- Only feeding when the CM300 run Material build up will damage the instrument's parts!
- Feed size must not exceed the maximum size of equipment Will make the CM300 stop functioning and may damage the device unit
- Use a slow and continuous feed . Excess feed will make the CM300 to stop running and will damage the machine parts
- During grinding, the dust may be run out from the feed hopper. When confronted with damage to the body of toxic or abrasive materials, use ventilation equipment. The risk of inhaling harmful dust!
- Many of the chemical materials are explosive mixture s by composition of the air after grinding . Please test your material specifications before grinding of the materials. Explosive hazard!



10.3 Replacement static knives


- Power cut off the instrument
- Open the door. Clean the grinding Chamber.
- Unscrew screw O and remove
- Replacement static knife N (4x)
- Tighten the screw O.

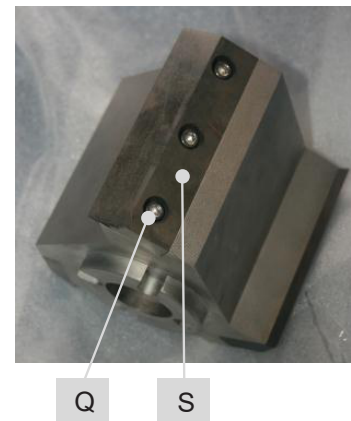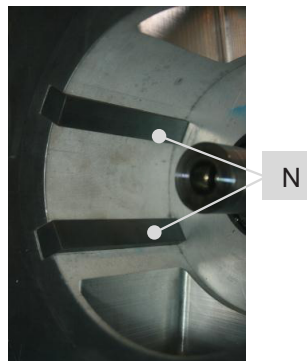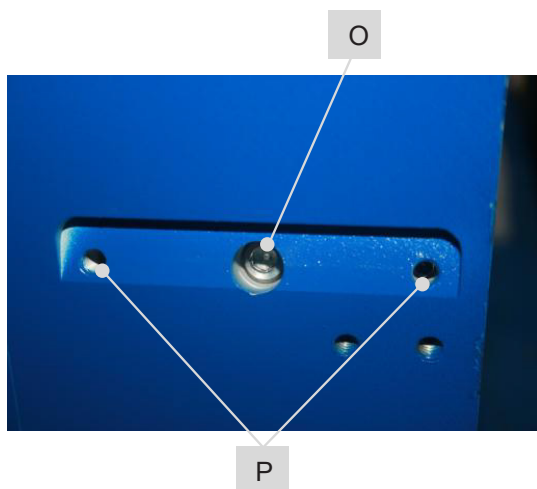

- Power cut off the instrument
- Open the door. Clean the grinding chamber.
- Pull out the rotor with the knives. Be care of personnel injuries by sharp knives!
- Unscrew screw Q (3x) and remove
- Replacement static knife N (3x)
- Tighten the screw O.


10.4 Adjust the gap

- Use a feeler gauge for detection gap S between the rotor knife and static knife N.
- Gap settings => 0.2mm, < 0.5mm.
- Reduce the gap:
 - Release screw O, screw P (2x) clockwise rotation, tighten the screw O.
 - Constantly check with a feeler until the gap to meet the requirements
 - Tighten screw O, applied torque = 7 N / m.
- Increasing the gap:
 - Release screw O,, symmetry parallel to loosen the screw P, counterclockwise rotation, tighten the screw O,
 - Continuously check with a feeler until the gap to meet the requirements.
 - Tighten screw O, applied torque = 7 N / m.

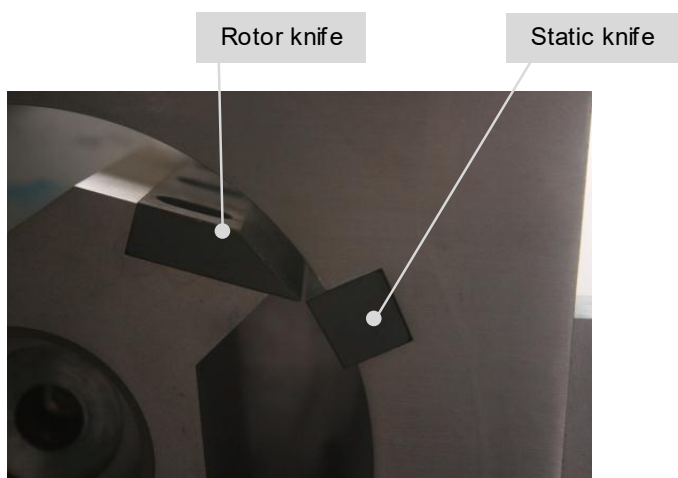

10.5 Adjusting by setting the rotor knife

The procedure as in **10.4** but then with the rotor knives. The hardened steel do have oval holes to change the position of the knife on the rotor block.

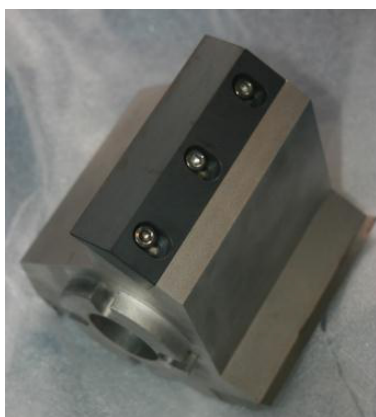


11 Preparation before grinding operation

- Disconnect the equipment power supply
- Insert the bottom sieve you selected
- Placed cutter
- Turn off the door grinding Chamber
- Fasten the collecting bucket

**11.1 Standard hopper**

The sample size hopper mainly for standard samples including:

- Cable waste
- Scrap carpets
- Straw, grass and similar materials
- Twigs and similar materials
- The material into the hopper
- Promote the use of plug material feed

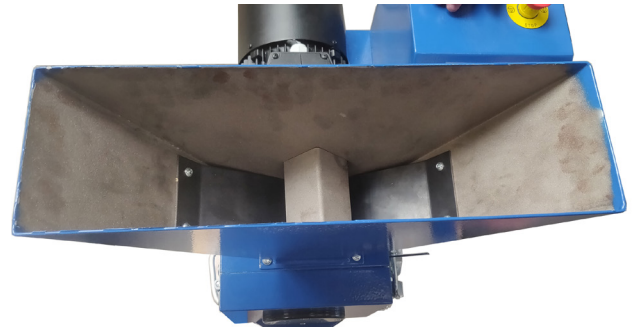**11.2 Start the instrument**

- Close the door of the grinding Chamber
- Set the grinding speed
- Pressing START key to start the equipment

11.3 Stop the instrument

- Pressing STOP key to stop the instrument

Handle long material, please bring protective gloves Or it may hurt the hand.

When push materials, only use the feed plug.

Use of other auxiliary equipment may result in damaged equipment accessories



12 Troubleshoot and Maintenance

This section helps identify and correct common issues during the operation. Always ensure the machine is switched off and disconnected from power before inspecting or performing maintenance.

12.1 Troubleshoot

Problem	Causes	Actions
Mill will not start (no display / no response)	No mains power; emergency stop engaged; door/cover not closed; blown fuse; power switch off	Check mains cable/outlet; reset emergency stop; close and latch door/cover; check/replace fuse per manual; switch on main power
Mill starts but stops immediately	Safety interlock open; overload protection triggered; motor thermal protection active	Ensure cover and collection unit are correctly installed; reduce load; let motor cool; restart after reset period
Motor overload / frequent tripping	Feeding too fast; material too hard/tough; wrong rotor speed; blunt knives; sieve blocked	Reduce feed rate; pre-cut/pre-cool sample; select appropriate speed; sharpen/replace knives; clean sieve and grinding chamber
Poor throughput / slow milling	Sieve aperture too small; blunt knives; wet/sticky sample; buildup in chamber; feed funnel bridging	Use larger sieve; sharpen/replace knives; pre-dry/pre-cool; clean chamber; use anti-bridging feed technique or smaller feed portions
Particle size too coarse	Sieve too large; worn/blunt knives; sample fed too quickly; rotor speed too low	Install smaller sieve; sharpen/replace knives; reduce feed rate; increase speed (within permitted range)
Excessive fines / dust	Sieve too small; rotor speed too high; brittle sample; over-milling	Use larger sieve; reduce speed; shorten milling time; use stepwise milling (coarse → fine)
Sample heats up / melts / smears	Speed too high; feed too fast; friction from blunt knives; thermoplastic sample; clogged sieve	Reduce speed/feed; sharpen knives; pre-cool sample (ice/dry ice/LN ₂ if allowed); pause to cool; clean/clear sieve
Sample sticks to chamber / forms lumps	Moist/oily sample; temperature rise; insufficient airflow/aspiration; incorrect sieve	Pre-dry/pre-cool; reduce speed; clean chamber; use suitable sieve; use cyclone/aspiration if available
Sieve clogs quickly	Wet/fibrous material; sieve aperture too small; too much feed; resinous material	Pre-dry/pre-cool; choose larger aperture; reduce feed; clean sieve more frequently; consider cutting in stages
Product backs up into feed hopper	Collection container not seated; filter/cyclone blocked; sieve clogged; airflow restricted	Fit collection unit correctly; empty/clean filter/cyclone; clean sieve; check seals/hoses if aspiration used
Excessive vibration	Rotor out of balance (buildup); damaged rotor; knives mounted incorrectly; loose fasteners; mill not level	Stop, clean rotor; inspect for damage; reinstall knives correctly; tighten fasteners to torque spec; level the mill
Unusual noise (scraping/metallic)	Knife contact due to wrong installation; foreign object; loose knife screws; damaged bearings	Stop immediately; inspect chamber; remove foreign object; reinstall/replace knives; tighten screws; service bearings if needed
Poor repeatability between runs	Inconsistent feed rate; mixed sample moisture; dull knives; different sieve; residue contamination	Standardize feed protocol; condition sample (drying/cooling); maintain knife sharpness; use same sieve; clean between runs
Cross-contamination between samples	Inadequate cleaning; residue trapped behind sieve/knives; shared collection container	Follow validated cleaning SOP; remove sieve and clean thoroughly; clean behind knives; use dedicated liners/containers
Chamber cannot be opened / door won't latch	Interlock misaligned; debris on seal; collection unit not seated; latch mechanism dirty	Remove debris; seat collection unit properly; clean latch; align interlock actuator; do not bypass interlocks
Rotor won't turn by hand (power off)	Jammed sample; foreign object; seized bearing; knives contacting	Disconnect power; clear jam; inspect for foreign objects; check knife clearance; service bearings if needed
Excessive wear of knives/sieves	Abrasive sample (minerals, glass); wrong knife material; running too fast; contaminant (sand/metal)	Use wear-resistant knives (e.g., hardened steel/tungsten carbide if offered); reduce speed; pre-screen contaminants; plan preventive replacements
Dust escaping / leakage	Seals worn; collection unit not sealed; aspiration not connected; filter overloaded	Check/replace seals; reseal and lock collection unit; connect aspiration/cyclone; clean/replace filter
Sample not cutting (fibers wrap around rotor)	Fibrous/tough material; wrong speed; dull knives; sieve too small; too large feed pieces	Pre-cut to shorter lengths; use suitable speed (often lower for tough fibers); sharpen knives; use larger sieve; feed smaller pieces
Rotor stalls with soft rubber/plastics	Material smears; too warm; speed too high; blunt knives	Pre-cool; reduce speed; sharpen knives; consider cryogenic milling if applicable

12.2 Maintenance and lubrication

Please view the service menu page for more information (chapter 9.6) Every 120 hours, the machine must be inspected for any form of wear and tear. Contact Laarmann for advice or replacement of wear parts or defective components.

General maintenance tasks

Daily	- Check the machine for visible damage or loose parts. - Ensure safety cover and interlocks are working properly. - Clean cutting chamber and surrounding area after each use.
Weekly	- Inspect power cable and plug for wear or damage. - Check that the drum rotates smoothly and evenly. - Verify that grinding media are in good condition (no deformation or cracks).
Monthly	- Check all bolts and fasteners for tightness. - Lubricate bearings or moving parts if applicable
Annually	- Have a qualified technician perform a full inspection of the motor, electrical components, and safety systems

12.3 Wear

The rotor and knives can become worn out, depending on the frequency of the milling operation and the milled material.

Always use original spare parts or replacements approved by the manufacturer. Using non-approved components may cause malfunction or void the warranty.

12.4 Cleaning

- Stop the CM300
- Turn off the power switch the back of the instrument control cabinet
- Open the door of the grinding chamber
- You can use brushes, industrial vacuum cleaners, compressed air for cleaning.
- If necessary, pull out by hand will be turn the cutter easy to clean.
- Clean the instrument only with a damp cloth.



Before cleaning the machine, unplug the mains cord from wall socket.
Use only a soft cloth moistened with water. Cleaning agents and solvents should not be used.



Do not clean the CM300 under running water. Danger to life through electric shock.



When the door is open, if the safety switch can not stop running the instrument, the instrument can not continue to use. A risk of personal injury!